



RESOURCE AND PATIENT MANAGEMENT SYSTEM

ICD-10 Update

(AUM)

Installation Guide and Release Notes

Version 26.0
October 2025

Office of Information Technology
Division of Information Technology

Table of Contents

1.0	Release Notes	1
1.1	Description of Changes	1
2.0	Installation Notes.....	5
2.1	Contents of Distribution	5
2.2	Required Resources	5
3.0	Installation Instructions	6
4.0	Sample Installation.....	7
4.1	Sample Cache Install.....	8
	Acronym List	13
	Contact Information	14

Preface

This manual contains the installation guide for the Indian Health Service (IHS) International Classification of Diseases, 10th Revision, Clinical Modification and Procedure Coding System (ICD-10-CM/PCS), v26.0.

1.0 Release Notes

This package will load new, modified, and invalid ICD-10 code updates for fiscal year (FY) 2026 for both diagnosis and operation/procedures codes. Additionally, the major diagnostic category and diagnostic-related grouper codes have been updated to fiscal year 2026. This update will ensure users have access to current codes for accurate and compliant visit information, claims reimbursement, and reporting.

1.1 Description of Changes

The 2026 ICD-10-CM files below contain information on the ICD-10-CM updates for FY 2026. These 2026 ICD-10-CM codes are to be used for discharges occurring from October 1, 2025 through September 30, 2026 and for patient encounters occurring from October 1, 2025 through September 30, 2026.

The 2026 ICD-10 Procedure Coding System (ICD-10-PCS) files below contain information on the ICD-10-PCS updates for FY 2026. These 2026 ICD-10-PCS codes are to be used for discharges occurring from October 1, 2025 through September 30, 2026.

The 2026 Medicare Severity Diagnosis-Related Groups (MS-DRG) files below contain information on the ICD-10 MS-DRG Grouper Logic Version 43.0 and Medicare Code Edits Version 43.0 updates and code lists effective for FY 2026.

As of AUM v18.0, the POA EXEMPT field (#1.9) is updated with data from the Centers for Medicare & Medicaid Services (CMS) for future integration.

As of FY 2019, CMS has removed the designation of codes when the primary diagnosis is its own Complications and Comorbidities (CC) or Major Complications and Comorbidities (MCC).

As of FY 2020, CMS has updated the Maternity Age range in the Medicare Code Edits v37 to be 9–64 years inclusive. Previously this range was 12–55 years inclusive.

As of FY 2025, CMS has deactivated the Sex conflict in the Medicare Code Edits v42 that detected inconsistencies between a patient's sex and any diagnosis or procedure on the patient's record.

- 2026 ICD-10-CM and PCS:
 - 2026 ICD-10-CM
(<https://www.cms.gov/medicare/coding-billing/icd-10-codes>)
 - 2026 ICD-10 PCS
(<https://www.cms.gov/medicare/coding-billing/icd-10-codes>)
 - ICD-10-CM/PCS MS-DRG v43.0

https://www.cms.gov/icd10m/FY2026-fr-v43-fullcode-cms/fullcode_cms/P0001.html

- Definitions of Medicare Code Edits v43.0
- ICD DIAGNOSIS File (#80):
 - .01 CODE NUMBER
 - 1.1 CODING SYSTEM
 - 1.3 UNACCEPTABLE AS PRINCIPAL DX
 - 1.9 POA EXEMPT
 - 10 SEX (Multiple-80.04)
 - .01 SEX EFFECTIVE DATE
 - 1 SEX
 - 11 AGE LOW (Multiple-80.011)
 - .01 AGE LOW EFFECTIVE DATE
 - 1 AGE LOW
 - 12 AGE HIGH (Multiple-80.012)
 - .01 AGE HIGH EFFECTIVE DATE
 - 1 AGE HIGH
 - 66 STATUS (Multiple-80.066)
 - .01 STATUS EFFECTIVE DATE
 - .02 STATUS
 - 67 DIAGNOSIS (Multiple-80.067)
 - .01 DIAGNOSIS EFFECTIVE DATE
 - 1 DIAGNOSIS
 - 68 DESCRIPTION (Multiple-80.068)
 - .01 DESCRIPTION EFFECTIVE DATE
 - 1 DESCRIPTION
 - 71 DRG GROUPER (Multiple-80.071)
 - .01 DRG GROUPER EFFECTIVE DATE
 - 1 DRG (Multiple-80.711)
 - .01 DRG
 - 72 MAJOR DIAGNOSTIC CATEGORY (Multiple-80.072)
 - .01 MDC EFFECTIVE DATE
 - 1 MDC

- 103 COMPLICATION/COMORBIDITY (Multiple-80.0103)
 - .01 CC EFFECTIVE DATE
 - 1 COMPLICATION/COMORBIDITY
- ICD OPERATION/PROCEDURE File (#80.1):
 - .01 CODE NUMBER
 - 1.1 CODING SYSTEM
 - 10 SEX (Multiple-80.11)
 - .01 SEX EFFECTIVE DATE
 - 1 SEX
 - 66 STATUS (Multiple-80.166)
 - .01 STATUS EFFECTIVE DATE
 - .02 STATUS
 - 67 OPERATION/PROCEDURE (Multiple-80.167)
 - .01 OPER/PROCEDURE EFFECTIVE DATE
 - 1 OPERATION/PROCEDURE
 - 68 DESCRIPTION (Multiple-80.168)
 - .01 DESCRIPTION EFFECTIVE DATE
 - 1 DESCRIPTION
 - 71 DRG GROUPE (Multiple-80.171)
 - .01 DRG GROUPE EFFECTIVE DATE
 - 1 MAJOR DIAGNOSTIC CATEGORIES (Multiple-80.1711)
 - .01 MAJOR DIAGNOSTIC CATEGORY
 - 1 DRG (Multiple-80.17111)
 - .01 DRG
- DRG File (#80.2):
 - .01 NAME
 - .06 SURGERY
 - 1 DESCRIPTION (Multiple-80.21)
 - .01 DESCRIPTION
 - 5 MDC#
 - 7 WEIGHT(nonAffil)
 - 10 AVG LENGTH OF STAY(days)
 - 14 ACTIVATION DATE

- 15 INACTIVE
- 16 INACTIVATION DATE
- 20 FISCAL YEAR WEIGHTS&TRIMS (Multiple-80.22)
 - 01 FISCAL YEAR WEIGHTS&TRIMS
 - 2.5 WEIGHT(nonAffil)
 - 4.5 AVG LENGTH OF STAY(days)
- 66 EFFECTIVE DATE (Multiple-80.266)
 - .01 EFFECTIVE DATE
 - .03 STATUS
 - .05 MDC#
 - .06 SURGERY
- 68 DESCRIPTION (VERSIONED) (Multiple-80.268)
 - .01 EFFECTIVE DATE
 - 1 DESCRIPTION (Multiple-80.2681)
 - .01 DESCRIPTION
 - 999999901 HCFA/IHS WEIGHT
 - 999999902 HCFA/IHS ARITH MEAN LOS
 - 999999903 HCFA/IHS GEOM MEAN
- MAJOR DIAGNOSTIC CATEGORY File (#80.3):
 - .01 NAME

2.0 Installation Notes

Prefix: AUM

Current version: 26.0

2.1 Contents of Distribution

Table 2-1: Contents of distribution—files and descriptions

File	Description
aum_2600.k	KIDS file to run the loading of the ICD-10 code updates
aum_260i.pdf	Installation Guide and Release Notes (this document)

2.2 Required Resources

- XU v8.0 p1018 or higher
- DI v22.0 p1018 or higher
- AICD v4.0
- AUM v25.0 or higher
- AUM v25.0 p2

3.0 Installation Instructions

In all User Class Identifications (UCIs) requiring ICD 2026 updates and in all databases running at least AUM v25.0 and AUM v25.0 p2:

1. Use the utility **D ^%GOGEN** to save the following Globals:
 - ^ICD9–ICD Diagnosis codes
 - ^ICD0–ICD Operation/Procedure codes
 - ^ICD–DRG codes
 - ^ICM–MDC codes
2. This release should be installed during off-peak hours.
3. Load the Distribution into KIDS using option **1** on the **Installation** menu. The distribution was released in a file named **aum_2600.k**.
4. Verify the load using option **2** on the **Installation** menu.
5. Consider using options **3** and **4** to print and compare the Transport.

<p>Note: For the next step, invoke the Kernel Installation & Distribution System menu by entering D ^XPDKRN from _PROGRAMMER MODE_.</p>
--

6. Install the distribution using option **6** on the **Installation** menu. Answer **NO** to the “Inhibit Logons” and “Disabling Scheduled Options” prompts as this build may have a longer-than-usual installation run time.
7. For results, use the KIDS **Install File Print** option to view and print the AUM 26.0 entry.

4.0 Sample Installation

The following subsection shows a sample installation. Screens may not be exactly the same during installation.

4.1 Sample Cache Install

```
D ^XUP

Setting up programmer environment
This is a PRODUCTION account.

Terminal Type set to: C-VT100

Select OPTION NAME: XPD MAIN          Kernel Installation & Distribution System

    Edits and Distribution ...
    Utilities ...
    Installation ...
    Patch Monitor Main Menu ...

Select Kernel Installation & Distribution System Option: INstallation

    1    Load a Distribution
    2    Verify Checksums in Transport Global
    3    Print Transport Global
    4    Compare Transport Global to Current System
    5    Backup a Transport Global
    6    Install Package(s)
        Restart Install of Package(s)
        Unload a Distribution

Select Installation Option: 1  Load a Distribution
Enter a Host File: C:\Patches\AUM 26.0\aum_2600.k

KIDS Distribution saved on Sep 11, 2025@11:49:48
Comment: ICD10 FY2026

This Distribution contains Transport Globals for the following Package(s):
    AUM 26.0
Distribution OK!

Want to Continue with Load? YES//
Loading Distribution...

Build AUM 26.0 has an Environmental Check Routine
Want to RUN the Environment Check Routine? YES// YES
    AUM 26.0
Will first run the Environment Check Routine, AUMI10

        Hello, SYSTEM MANAGER

        Checking Environment for ICD UPDATE V 26.0.

        Need at least XU v8.0 p1018.....XU v8.0 p1020 Present
```

```

Need at least DI v22.0 p1018.....DI v22.0 p1021 Present
Need at least AICD v4.0.....AICD v4.0 Present
Need at least AUM v25.0.....AUM v25.0 Present
Requires AUM v25.0 p2.....Present

ENVIRONMENT OK.
Use INSTALL NAME: AUM 26.0 to install this Distribution.

1      Load a Distribution
2      Verify Checksums in Transport Global
3      Print Transport Global
4      Compare Transport Global to Current System
5      Backup a Transport Global
6      Install Package(s)
      Restart Install of Package(s)
      Unload a Distribution

Select Installation Option:

Edits and Distribution ...
Utilities ...
Installation ...
Patch Monitor Main Menu ...

Select Kernel Installation & Distribution System Option:

Do you really want to halt? YES//

Logged out at Sep 12, 2025 11:30 am
GOLDB>D ^XPDKRN

KIDS 8.0

Select KIDS OPTION: INSTALL PACKAGE(S)
Select INSTALL NAME: AUM 26.0          Loaded from Distribution    9/12/25@11:30:44
=> ICD10 FY2026 ;Created on Sep 11, 2025@11:49:48

This Distribution was loaded on Sep 12, 2025@11:30:44 with header of
ICD10 FY2026 ;Created on Sep 11, 2025@11:49:48
It consisted of the following Install(s):
AUM 26.0
Checking Install for Package AUM 26.0
Will first run the Environment Check Routine, AUMI10

Hello, SYSTEM MANAGER

Checking Environment for ICD UPDATE V 26.0.

```

```

Need at least XU v8.0 p1018.....XU v8.0 p1020 Present
Need at least DI v22.0 p1018.....DI v22.0 p1021 Present
Need at least AICD v4.0.....AICD v4.0 Present
Need at least AUM v25.0.....AUM v25.0 Present
Requires AUM v25.0 p2.....Present

ENVIRONMENT OK.

Install Questions for AUM 26.0

Incoming Files:

9003103 AUM DATA (including data)
Note: You already have the 'AUM DATA' File.
I will OVERWRITE your data with mine.

Want KIDS to INHIBIT LOGONS during the install? NO//

Install Started for AUM 26.0 :
Sep 12, 2025@11:30:45

Build Distribution Date: Sep 11, 2025

Installing Routines:
Sep 12, 2025@11:30:45

Running Pre-Install Routine: PRE^AUMI10

Installing Data Dictionaries:
Sep 12, 2025@11:30:46

Installing Data:
Sep 12, 2025@11:31:18

Running Post-Install Routine: POST^AUMI10
-----

MAJOR DIAGNOSTIC CATEGORY FILE
MDC NAME
=== =====

UNMODIFIED.....27
TOTAL.....27
-----

DRG FILE
DRG DESCRIPTION MDC SURG WEIGHT LOS
=== =====
MODIFIED.....769
DESCRIPTION..11

```

```

LOS.....500
MDC.....5
STATUS.....5
SURGERY.....2
WEIGHT.....768
UNMODIFIED.....3
TOTAL.....772
-----

ICD10 DIAGNOSIS FILE
      CODE      STA SHORT LONG SEX MDC CC  LOW HIGH DRG PDX POA
      ====      ==  ===== ===== ==  ==  ==  ==  =====
NEW.....487
MODIFIED.....218
  CC.....2
  DRG.....171
  LONG.....38
  MDC.....4
  PDX.....7
  SHORT.....38
INACTIVE.....28
UNMODIFIED.....97465
TOTAL.....98198
-----

ICD10 OPERATION/PROCEDURE FILE
      CODE      STA SHORT LONG SEX MS-DRG
      ====      ==  ===== ===== ==  =====
NEW.....156
MODIFIED.....367
  DRG.....367
INACTIVE.....27
UNMODIFIED.....79518
TOTAL.....80068
-----

ICD-10 Table Update Summary
                                New      Mod      Ina      Error      Unmod      Total
MDC                                0        0        0        0        27        27
DRG                                0       769        0        0         3       772
ICD-10-CM                        487       218       28        0      97465     98198
ICD-10-PCS                       156       367       27        0      79518     80068
-----


Updating Routine file...

Updating KIDS files...

AUM 26.0 Installed.
          Sep 12, 2025@11:32:37

NO Install Message sent

```

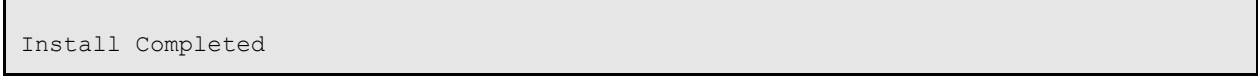

Install Completed

Figure 4-1: Sample install

Acronym List

Acronym	Term Meaning
CC	Complications and Comorbidities
CMS	Centers for Medicare & Medicaid Services
FY	Fiscal Year
ICD-10-CM/PCS	International Classification of Diseases, 10th Revision, Clinical Modification and Procedure Coding System
IHS	Indian Health Service
MCC	Major Complications and Comorbidities
MS-DRG	Medicare Severity Diagnosis-Related Groups
RPMS	Resource and Patient Management System
UCI	User Class Identification

Contact Information

If you have any questions or comments regarding this distribution, please contact the IHS IT Service Desk.

Phone: (888) 830-7280 (toll free)

Web: <https://www.ihs.gov/itsupport/>

Email: itsupport@ihs.gov